


Synergy effect

RM | Microneedle

Needles are adjustable from depth 0.5mm to 3.5mm, and provide insulated and non-insulated needle effect.

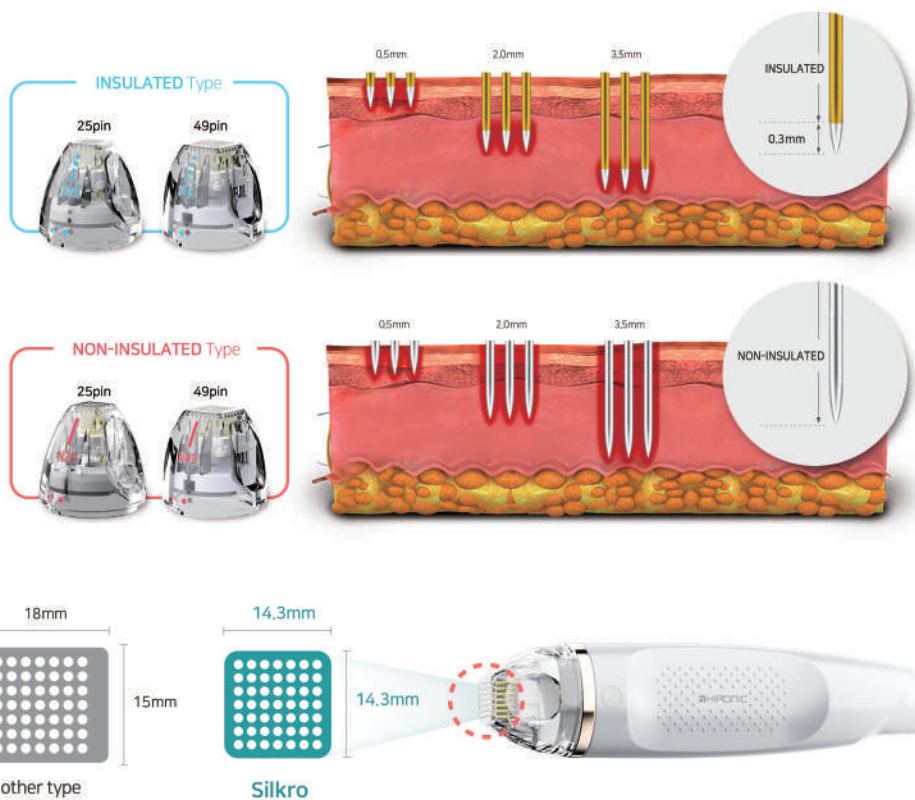

Wide viewing angle

With 25% smaller surface area of tip size, it can minimize overlap during the treatment and provide a wide viewing angle for more accurate and precise procedures.

PS | PLASMA

Plasma is effective in whitening, sterilization, tightening, and helps absorb the best cosmetic ingredient 466% more. Its unique features and technology are proven through clinical trials.



Silkro Specification

Handpiece	Radio Frequency	RF - Micro Needle	25pin	Insulated
			49pin	Non-insulated
		RF - Circle*	30V ~ 360V	
		RF - Needle 1pin*	10V ~ 168V	
		RF - Vaginal*	10V ~ 168V	
	PLASMA*	Level 1 ~ 5		
Shot count		2,000 shots		
Penetration Depth		0.5 ~ 3.5mm (0.1 step)		
Repetition		Single / 0.2 / 0.5 / 0.8 / 1.0 / 2.0 (sec)		
RF Intensity		Level 1 ~ 10		
Treating size		25PIN (8×8mm) / 49PIN (12×12mm)		
Frequency		2MHz		
RF Emission Time		50 ~ 950ms		
Convenience system		12.1" Wide LCD		
		5 language support (EN, CN, KR, ES, AR)		
		Voice Guidance		
		Selective mode		
		Memory function		
General		Air vent hole		
		Max Output	25W (25pin) / 50W (49pin)	
		Dimensions	654(W) x 768(L) x 1287(H)	
		Weight	43kg	

About Hironic?

HIRONIC Co., Ltd. is one of the best medical & aesthetic companies listed on the Korean stock market (KOSDAQ), and is leading the global beauty medical device market based on its management ideology as a "company that creates a happy world" with prime technology. HIRONIC continues to invest in R&D to deliver new value to our customers focusing on innovative and differentiated products. HIRONIC released innovative equipment in various beauty medical markets such as HIFU (MFU), RF, Laser, and Plasma devices. HIRONIC will continue to think about customers first and move towards future innovate technology and services.



19F U-Tower, 767, Sinsu-ro, Yongin-si, Gyeonggi-do, 16827 korea

Sales Support. Tel. +82-31-525-7000 Fax. +82-31-525-7010 E-mail. info@hironic.com www.hironic.com
Customer Service. Tel. +82-31-525-7600 E-mail. ocs@hironic.com



Hironic has achieved ISO 13485 and Korean GMP.

SK(CE)_02302_V001



For Export

RF Total Solution
silkro





Innovative Combination of RF and PLASMA


Silkro is the first innovative RF device unified with PLASMA technology

It can maximize the skin treatment effect from the epidermis to the dermis and delivers optimal result such as acne, wrinkle, tightening and skin rejuvenation overall.

RN | RF-Needle RF Monopolar 1pin needle

Selective sebaceous glands destruction

RF Needle selectively destroys sebaceous glands. Not damaging the epidermis, treating acne fundamentally, and preventing recurrence.

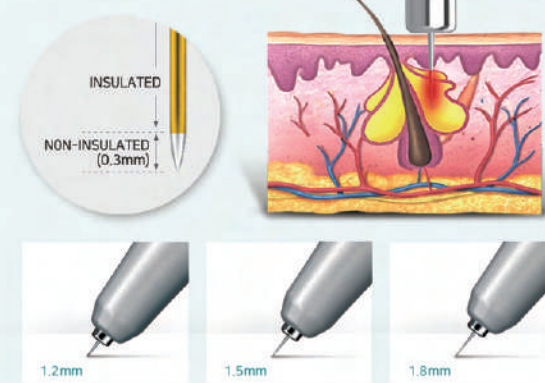


Thermal effect from 0.9mm, 1.2mm, 1.5mm

RC | RF-Circle RF Monopolar Circle

Bulk heating before and after treatment


With the small RF ceramic, transfer RF energy to the skin that is difficult to target by needle. Use for before and after acne treatment and general rejuvenation, and for women's vulvar treatment.

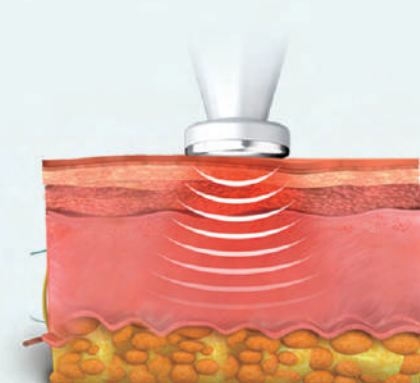

RV | RF-Vaginal RF Monopolar Bar type

Special Vaginal Treatment

Convenient radiofrequency treatment specially for women. Non-invasive, non-ablative, non-surgical, and no anesthesia is required. Specially ergonomic-designed handpiece for treatments like vaginal rejuvenation & tightening and colpoxerosis.

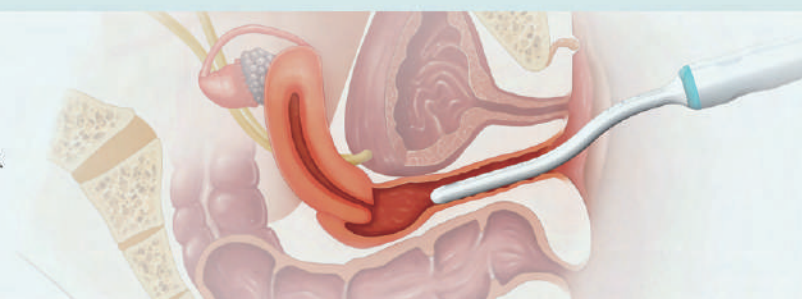

Treatment Lesion

RF Treatment

Scar treatment, Stretch marks, Skin lifting, Overall rejuvenation, Lifting, Acne, Black and white comedo, Pore, Fine wrinkle, Skin revitalization, Resurfacing, Hyperhidrosis, Tightening.

PLASMA Treatment

Skin Texture, Toning effect, Whitening, Dryness, Pore density, Hydration, Wrinkles, Drug delivery.


RF for Women

Vaginal Tightening, Vulva whitening, Vaginal rejuvenation.

Convenience

1 Smart Integrated System

Self-device monitoring & diagnosis through online and enabling real-time troubleshooting without visiting.

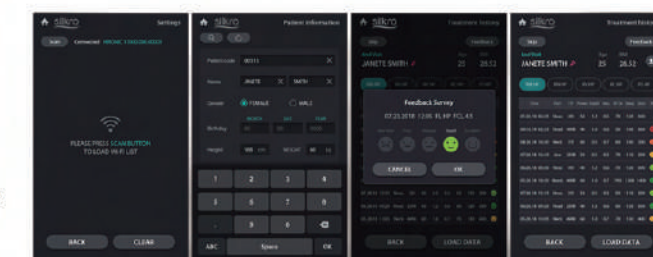

R.M.S (Remote Maintenance System)

Through accumulative treatment data, providing the best and safety protocol for each lesion for giving utmost guide to users.

- Consumable replacement notice
- Protocols saved
- Self diagnosis and inspection
- Malfunction prevention and assist functions
- Checking of accumulated shots and device use is available (Accumulated shots, used time and etc.)
- Remote update of SW/FW is available

Real-time Device monitoring

Patients clinical data

Pre-A/S



2 Big Data System

Through accumulative treatment data, providing the best and safest protocol for each lesion, and giving utmost guidance to users.



12.1" Wide LCD Max 75° Angle



Sensor Technology

Handpiece automatic recognition system



Back handle Convenient mobility



Air vent hole Fast and powerful air circulator



LED guide mode Red, Yellow, White, Green

